



MINING DEALER BEST PRACTICE

S.O.STM Lab Data Entry Improvement using New Electronic System

Maintenance and
Repair

Application

Component Life
Management

Component
Renewal (CRC)

MARC
Management

1.0	Introduction	1
2.0	Best Practice Description.....	1
3.0	Implementation Steps	2
4.0	Benefits.....	5
5.0	Resources Required	6
6.0	Supporting Attachments / References	6
7.0	Related Best Practices	6
8.0	Acknowledgements.....	6

DISCLAIMER: The information and potential benefits included in this document are based upon information provided by one or more Cat[®] dealers, and such dealer(s) opinion of “Best Practices.” Caterpillar makes no representation or warranty about the information contained in this document or the products referenced herein. Caterpillar welcomes additional “Best Practice” recommendations from our dealer network.

January 2015
1411-3.03-1357

1.0 Introduction

Oil Cleanliness Analysis is one of the primary functions of S.O.STM program. This analysis is carried out for 80% of all incoming samples to S.O.STM laboratory. This test determines the presence of particles in the oil, which are classified by many symbols.


Before, we had manually written down all results of this test. Manual Data entry increases the time for Oil Cleanliness Analysis. It is double work, when we first write the results on paper and then we transfer these same results to S.O.STM Service manager. Using the paper-based information makes it inconvenient to keep this information and search for these results.

This Best Practice is aimed to improve the procedures of oil cleanliness analysis with using an electronic device.

2.0 Best Practice Description

Borusan Makina Kazakhstan has developed the software application which improves the process of visual analysis of oils.

Before

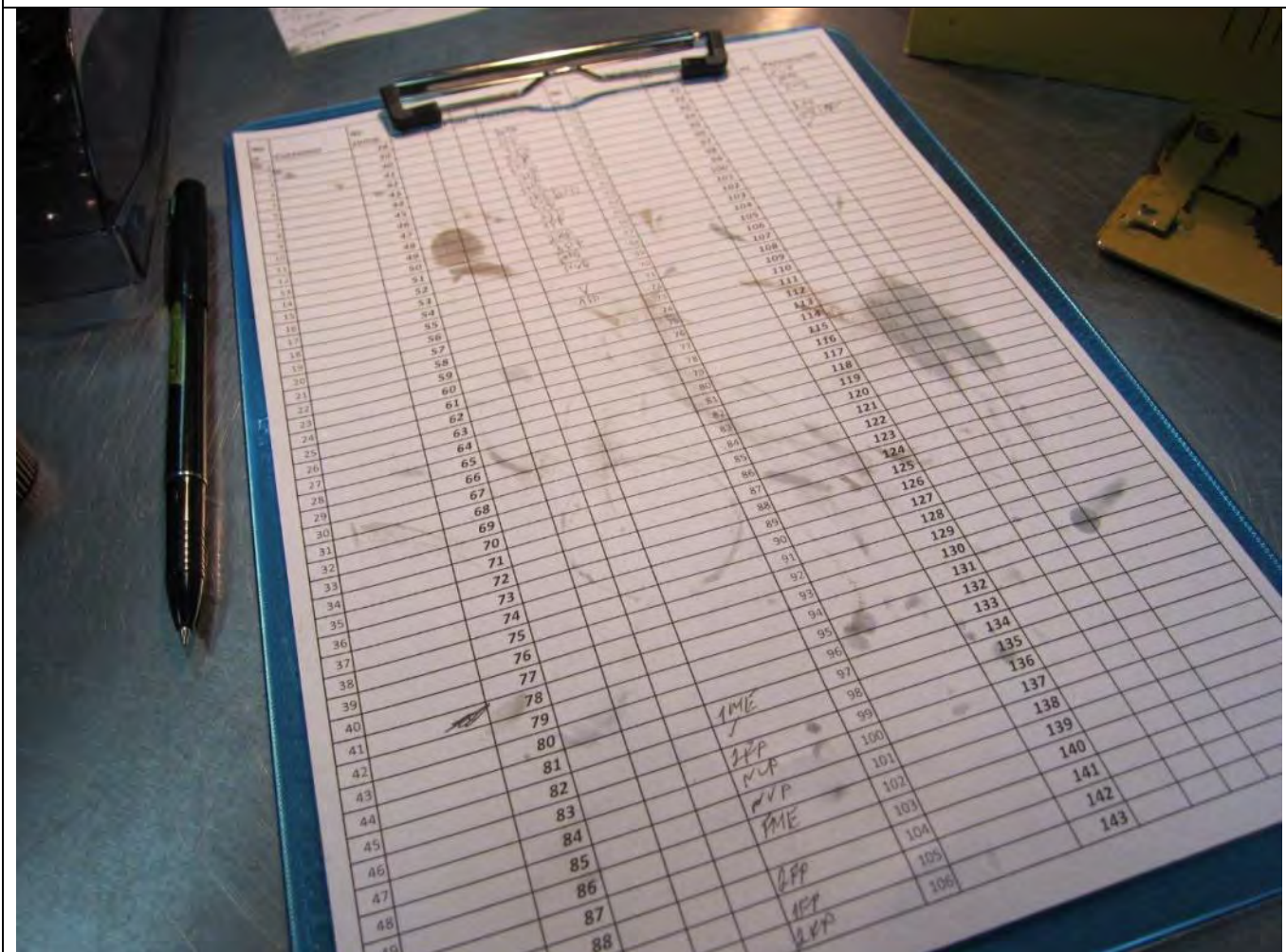

THE INFORMATION HEREIN MAY NOT BE COPIED OR TRANSMITTED TO OTHERS WITHOUT THE WRITTEN PERMISSION OF CATERPILLAR

SOS Lab Data Entry Improvement using New
Electronic System

DATE
2/16/2015

CHG
NO
00

NUMBER
1411-3.03-1357

After



3.0 Implementation Steps

The new process of analysis.

1. First, prior to the analysis, the employee should open the software on PC, which has been developed by IT specialist.


THE INFORMATION HEREIN MAY NOT BE COPIED OR TRANSMITTED TO OTHERS WITHOUT THE WRITTEN PERMISSION OF CATERPILLAR

SOS Lab Data Entry Improvement using New
Electronic System

DATE
2/16/2015

CHG
NO
00

NUMBER
1411-3.03-1357

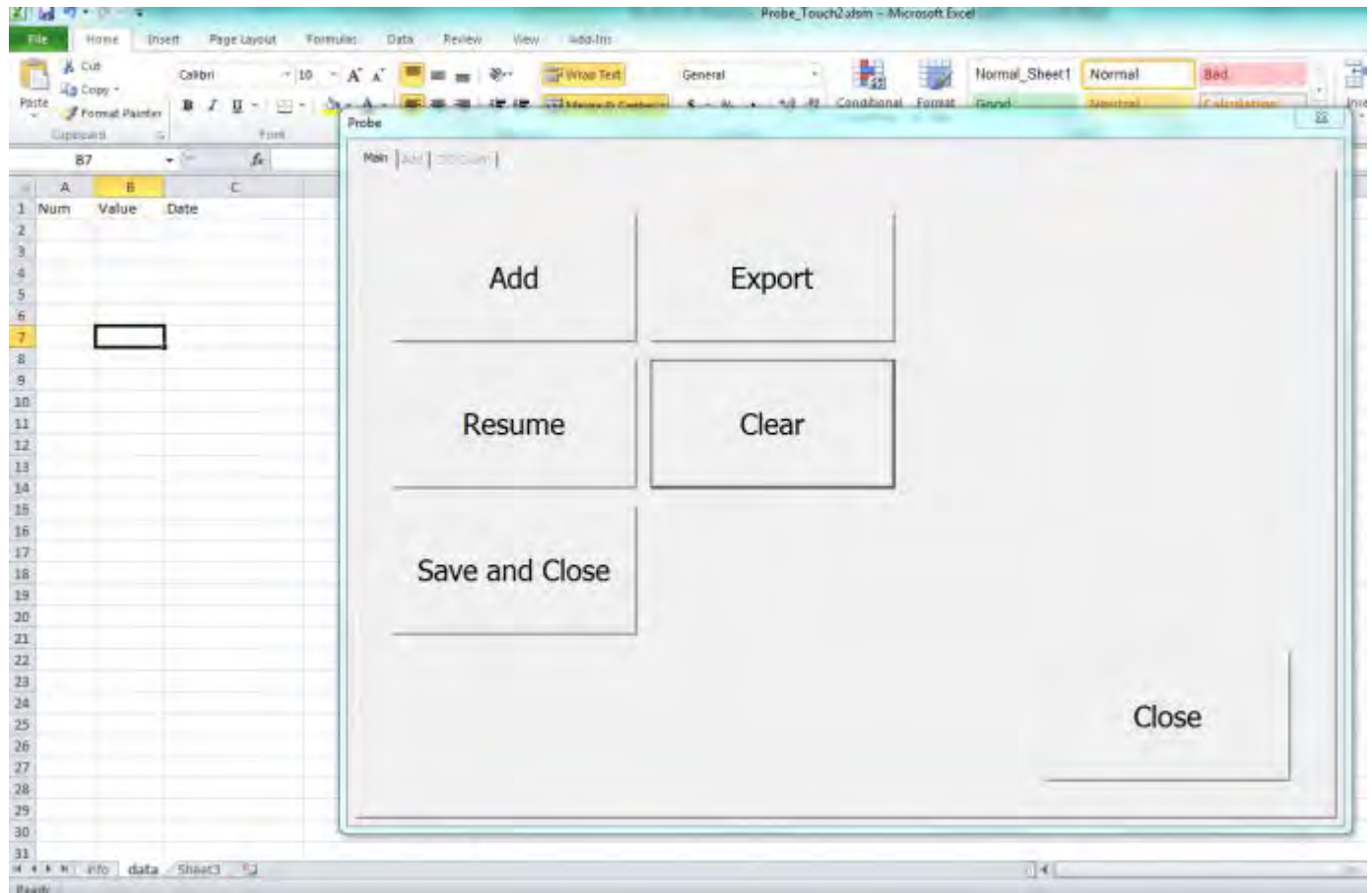

2. Then, using Tablet, we connect to the PC though the Remote Desktop Control Program.
3. Make visual evaluations of bottle cap.

THE INFORMATION HEREIN MAY NOT BE COPIED OR TRANSMITTED TO OTHERS WITHOUT THE WRITTEN PERMISSION OF CATERPILLAR

SOS Lab Data Entry Improvement using New
Electronic System


DATE
2/16/2015

CHG
NO
00

NUMBER
1411-3.03-1357



4. If the sample is clear, use the ISO button, and after completion of the analysis, the software prepares the list of samples, which should be analyzed on the Particle Counter.

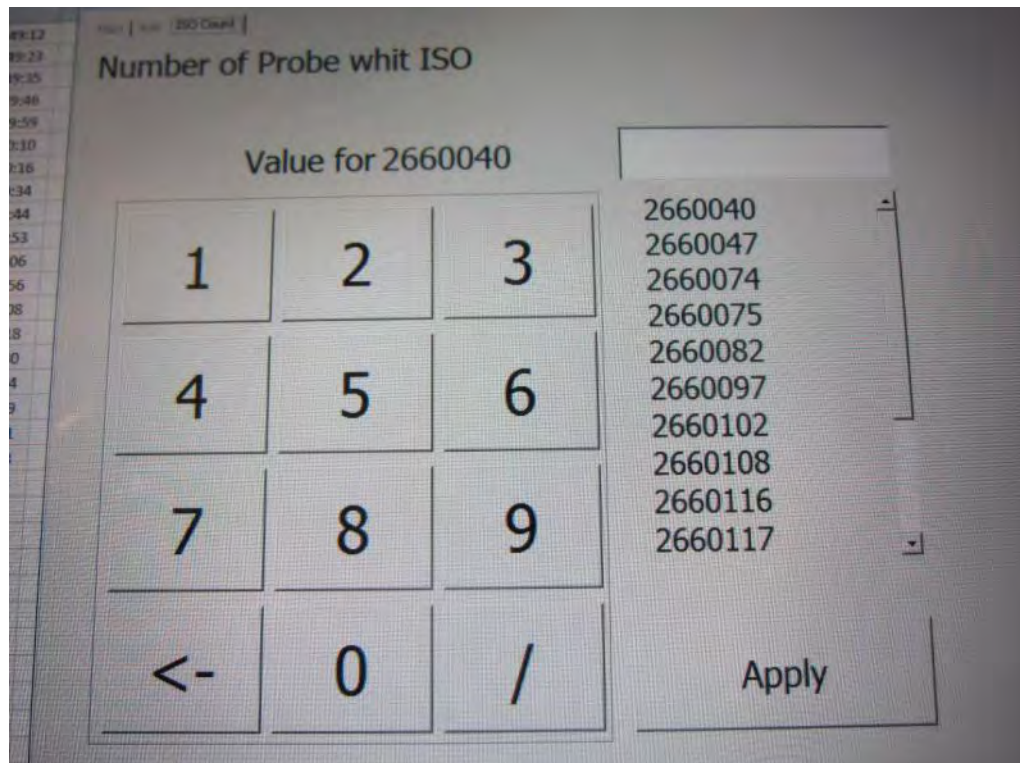

THE INFORMATION HEREIN MAY NOT BE COPIED OR TRANSMITTED TO OTHERS WITHOUT THE WRITTEN PERMISSION OF CATERPILLAR


SOS Lab Data Entry Improvement using New
Electronic System

DATE
2/16/2015

CHG
NO
00

NUMBER
1411-3.03-1357


- After full completion, the employee pushes 'Export' button, and program automatically create the txt file with necessary fields to upload them in Dataxfer folder.

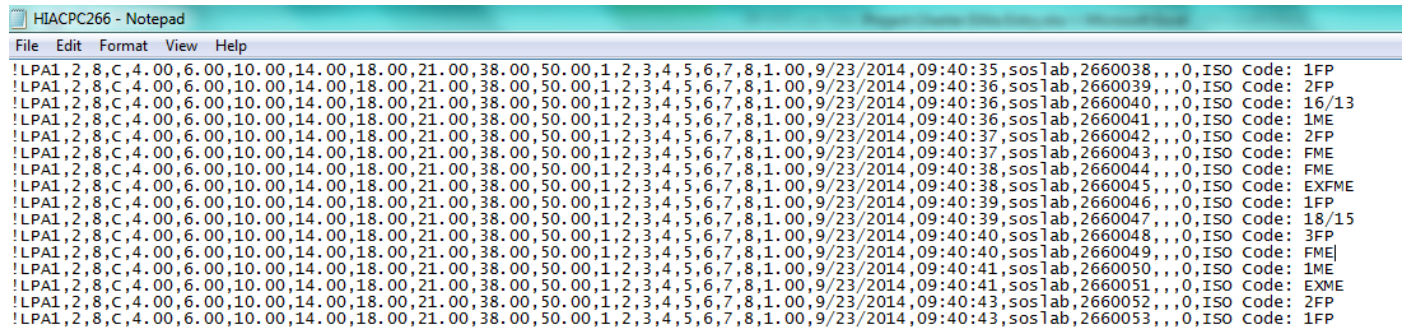

```

HIACPC266 - Notepad
File Edit Format View Help
!LPA1,2,8,C,4.00,6.00,10.00,14.00,18.00,21.00,38.00,50.00,1,2,3,4,5,6,7,8,1.00,9/23/2014,09:40:35,soslab,2660038,,,0,ISO Code: 1FP
!LPA1,2,8,C,4.00,6.00,10.00,14.00,18.00,21.00,38.00,50.00,1,2,3,4,5,6,7,8,1.00,9/23/2014,09:40:36,soslab,2660039,,,0,ISO Code: 2FP
!LPA1,2,8,C,4.00,6.00,10.00,14.00,18.00,21.00,38.00,50.00,1,2,3,4,5,6,7,8,1.00,9/23/2014,09:40:36,soslab,2660040,,,0,ISO Code: 16/13
!LPA1,2,8,C,4.00,6.00,10.00,14.00,18.00,21.00,38.00,50.00,1,2,3,4,5,6,7,8,1.00,9/23/2014,09:40:36,soslab,2660041,,,0,ISO Code: 1ME
!LPA1,2,8,C,4.00,6.00,10.00,14.00,18.00,21.00,38.00,50.00,1,2,3,4,5,6,7,8,1.00,9/23/2014,09:40:37,soslab,2660042,,,0,ISO Code: 2FP
!LPA1,2,8,C,4.00,6.00,10.00,14.00,18.00,21.00,38.00,50.00,1,2,3,4,5,6,7,8,1.00,9/23/2014,09:40:37,soslab,2660043,,,0,ISO Code: FME
!LPA1,2,8,C,4.00,6.00,10.00,14.00,18.00,21.00,38.00,50.00,1,2,3,4,5,6,7,8,1.00,9/23/2014,09:40:38,soslab,2660044,,,0,ISO Code: FME
!LPA1,2,8,C,4.00,6.00,10.00,14.00,18.00,21.00,38.00,50.00,1,2,3,4,5,6,7,8,1.00,9/23/2014,09:40:38,soslab,2660045,,,0,ISO Code: EXFME
!LPA1,2,8,C,4.00,6.00,10.00,14.00,18.00,21.00,38.00,50.00,1,2,3,4,5,6,7,8,1.00,9/23/2014,09:40:39,soslab,2660046,,,0,ISO Code: 1FP
!LPA1,2,8,C,4.00,6.00,10.00,14.00,18.00,21.00,38.00,50.00,1,2,3,4,5,6,7,8,1.00,9/23/2014,09:40:39,soslab,2660047,,,0,ISO Code: 18/15
!LPA1,2,8,C,4.00,6.00,10.00,14.00,18.00,21.00,38.00,50.00,1,2,3,4,5,6,7,8,1.00,9/23/2014,09:40:40,soslab,2660048,,,0,ISO Code: 3FP
!LPA1,2,8,C,4.00,6.00,10.00,14.00,18.00,21.00,38.00,50.00,1,2,3,4,5,6,7,8,1.00,9/23/2014,09:40:40,soslab,2660049,,,0,ISO Code: FME]
!LPA1,2,8,C,4.00,6.00,10.00,14.00,18.00,21.00,38.00,50.00,1,2,3,4,5,6,7,8,1.00,9/23/2014,09:40:41,soslab,2660050,,,0,ISO Code: 1ME
!LPA1,2,8,C,4.00,6.00,10.00,14.00,18.00,21.00,38.00,50.00,1,2,3,4,5,6,7,8,1.00,9/23/2014,09:40:41,soslab,2660051,,,0,ISO Code: EXME
!LPA1,2,8,C,4.00,6.00,10.00,14.00,18.00,21.00,38.00,50.00,1,2,3,4,5,6,7,8,1.00,9/23/2014,09:40:43,soslab,2660052,,,0,ISO Code: 2FP
!LPA1,2,8,C,4.00,6.00,10.00,14.00,18.00,21.00,38.00,50.00,1,2,3,4,5,6,7,8,1.00,9/23/2014,09:40:43,soslab,2660053,,,0,ISO Code: 1FP
  
```

- The Results are automatically uploaded to the appropriate samples in S.O.STM Service Manager Program.

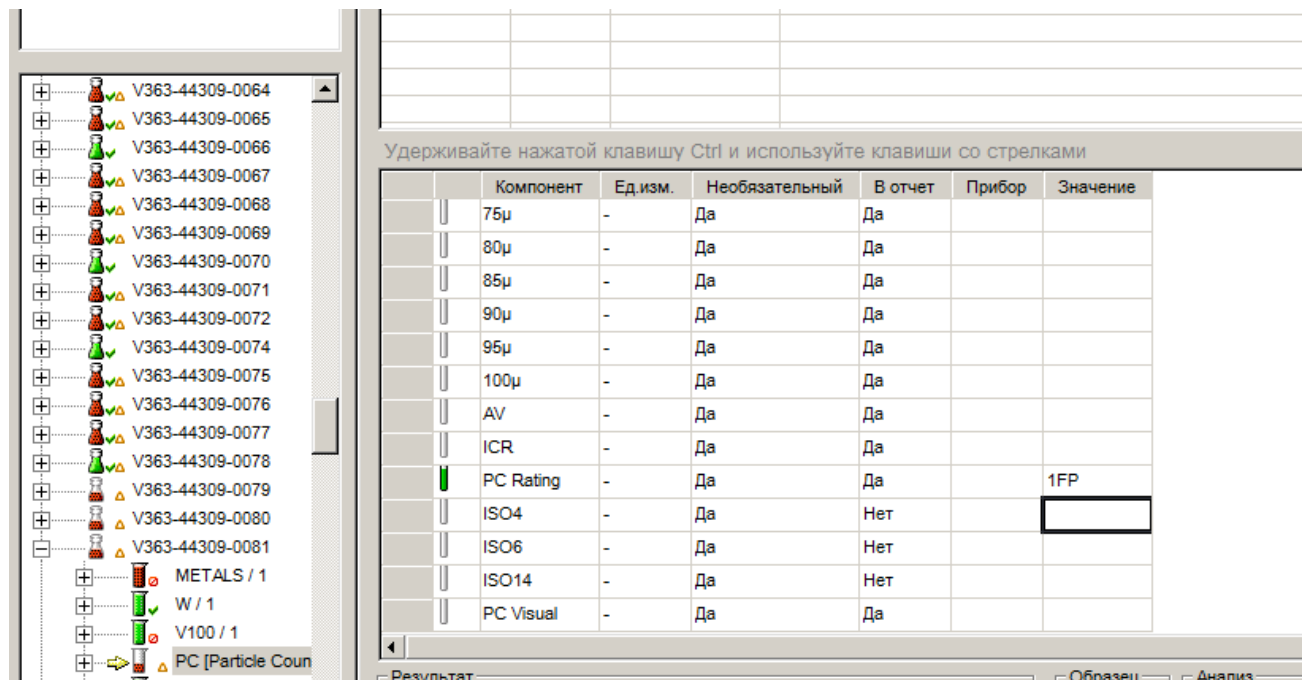

Компонент	Ед.изм.	Необязательный	В отчет	Прибор	Значение
75μ	-	Да	Да		
80μ	-	Да	Да		
85μ	-	Да	Да		
90μ	-	Да	Да		
95μ	-	Да	Да		
100μ	-	Да	Да		
AV	-	Да	Да		
ICR	-	Да	Да		
PC Rating	-	Да	Да		1FP
ISO4	-	Да	Нет		
ISO6	-	Да	Нет		
ISO14	-	Да	Нет		
PC Visual	-	Да	Да		

4.0 Benefits

- Data entry convenience
- Freedom of paper data storage
- Software generates the list of samples for Particle Counting
- Automatic export to the S.O.STM Service Manager program

THE INFORMATION HEREIN MAY NOT BE COPIED OR TRANSMITTED TO OTHERS WITHOUT THE WRITTEN PERMISSION OF CATERPILLAR			
SOS Lab Data Entry Improvement using New Electronic System	DATE 2/16/2015	CHG NO 00	NUMBER 1411-3.03-1357

5.0 Resources Required

- 1 Tablet PC
- Dealer must use LIMS software - S.O.STM Services Manager Program
- Specifically developed application

6.0 Supporting Attachments / References

If another dealer would like to replicate this software, they may contact Gulnaz Mansurova at Borusan for details.

Gulnaz Mansurova, CRC Process Engineer
 +7 (721) 248 26 26 ext. 117
 Oktyabrsky district, Karaganda-Temirtau 12 km,
 opposite "Ispat-Karmet" tennis court
 100005 Karaganda, Kazakhstan
gmansurova@borusan.com

7.0 Related Best Practices

N/A

8.0 Acknowledgements

The idea of creating & implementing this Best Practice belongs to:

Laura IVANSKAYA
 S.O.STM Lab Analyst
 Borusan Makina Kazakhstan LLP
livanskaya@borusan.com

Alexander PYAK
 IT Helpdesk Specialist
 Borusan Makina Kazakhstan LLP
apyak@borusan.com

Zhumat NURPEISSOV
 S.O.STM Lab Analyst
 Borusan Makina Kazakhstan LLP
znurpeissov@borusan.com

Natalya KOROLKOVA
 S.O.STM Lab Supervisor
 Borusan Makina Kazakhstan LLP
nkorolkova@borusan.com

This Best Practice was written by:

Natalya KOROLKOVA
 S.O.STM Lab Supervisor
 Borusan Makina Kazakhstan LLP
nkorolkova@borusan.com
 +7 (701) 763 97 10

THE INFORMATION HEREIN MAY NOT BE COPIED OR TRANSMITTED TO OTHERS WITHOUT THE WRITTEN PERMISSION OF CATERPILLAR

SOS Lab Data Entry Improvement using New Electronic System	DATE 2/16/2015	CHG NO 00	NUMBER 1411-3.03-1357
--	-------------------	-----------------	--------------------------